



FULL USER MANUAL

About GLO-ME

Your decision to purchase the **GLO-ME** Gentle Diamond Peel To Go has set you on the path to healthier and glowing skin. In order to fully benefit from your new **GLO-ME** device, it is important to carefully read this user manual. **GLO-ME** is a portable gentle diamond peel device which can simply and safely be used– anywhere, anytime.

CONTENTS	PAGE
1. Warnings & contraindications	1
2. Package contents	2
3. Directions for use	3-4
4. Cleaning & maintenance	5
5. Customer support & warranty	5-6
6. Technical specifications	6



WARNINGS

To reduce the risk of injury or product damage:

- Use this device for its intended purpose only, as described in this manual.
 - Keep out of the reach of children.
 - Do not store your GLO-ME in an area of high temperature or humidity. Prolonged exposure to heat or humidity may damage the device.
 - Do not submerge the unit in water or any other liquid. Hands must be completely dry when operating the device.
 - Sudden impact or electrical surges may damage the device. Do not use the device until it has been checked and/or serviced by a technician.
 - Do not attempt to open or disassemble the device, except for replacing treatment tips and filters. Do not use the device if disassembled.
 - Work only on intact skin. Do not use directly above or below your eyelid area, on the eye socket and on the lips, ears, breast and genital areas.
-


CONTRAINDICATIONS

Do not use the GLO-ME without consulting a physician if any of the following conditions apply to you:


- Use of **chemical peel** products, i.e. Alpha and Beta Hydroxy acid, Retin-A, 2 weeks before or after using GLO-ME or before skin is completely healed.
- Any injured active skin condition within treatment area such as sores, psoriasis, eczema, rash, undiagnosed lesions, active rosacea, warts, sunburn, inflamed acne or capillary area.
- Current or history of skin cancer or pre-malignant moles.
- Impaired immune system.
- Recent herpes outbreak.
- Unstable diabetes mellitus.
- Vascular lesions.
- Use of oral anti-coagulants (blood thinners) or Isotretinoin.
- If you are taking medications for skin conditions without your doctor's consent.
- Pregnancy.
- Excessive fresh skin tan (within a few days prior to treatment).
- If you are seeing a doctor for any skin condition.
- If you have very sensitive skin.

PACKAGE CONTENTS


- GLO-ME device (Fig 1)
- 2 types diamond peeling tips (Fig 2)
- 50 disposable filters (Fig 3)
- 2 batteries (1.5V, AAA size) (Fig 4)
- Specially designed carrying case (Fig 5)

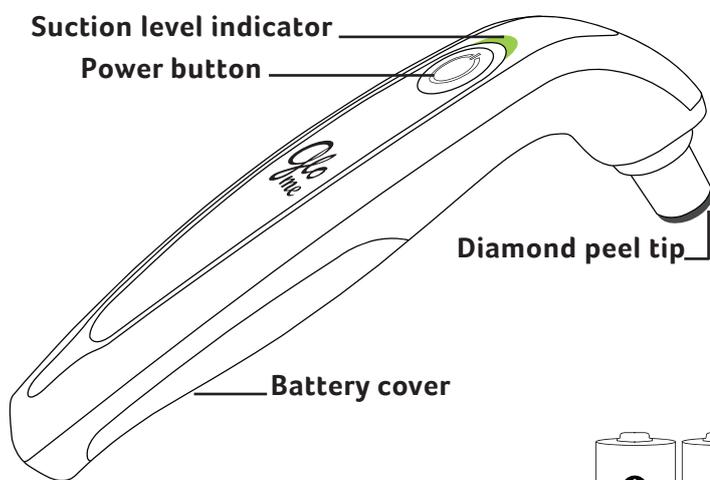


(Fig 1)

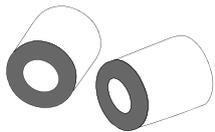


(Fig 2)

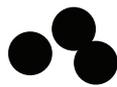

(Fig 3)

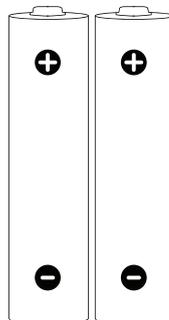

(Fig 4)

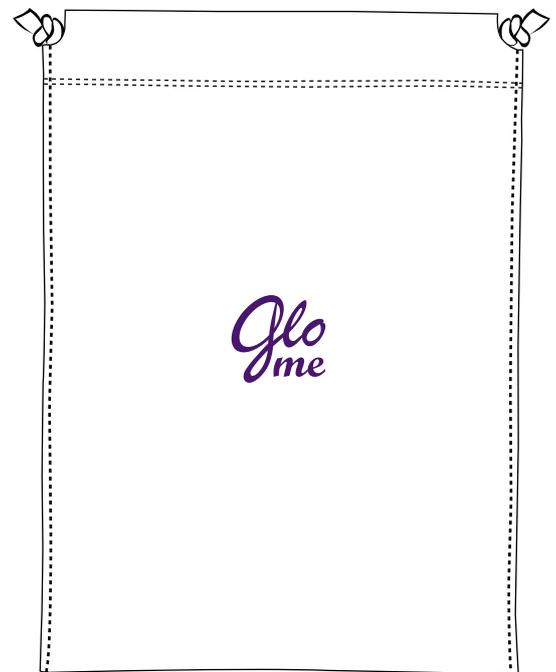

(Fig 5)


DIRECTIONS FOR USE

1. Open the battery cover and install two 1.5V AAA size batteries in the proper +/- direction. Then close the batteries cover and lock it.
2. Wash your skin to remove make-up, oils and dirt.
3. Make sure your skin is dry.
4. Place a new filter on the knob (Fig 6).
5. Select the type of tip* and connect it to the device by twisting and tightening it firmly against the device: (Fig 7).


Tip type*

White For **sensitive skin**

Silver For **normal skin**

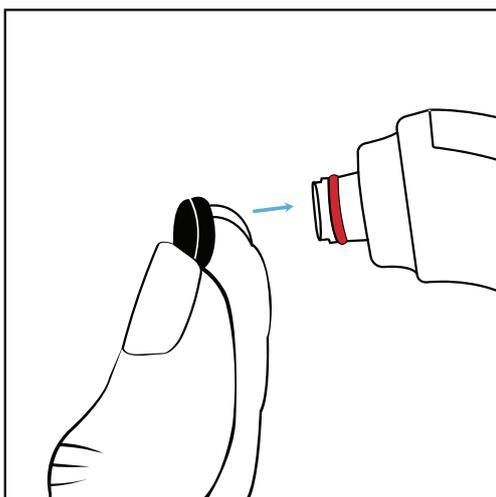


(Fig 6)

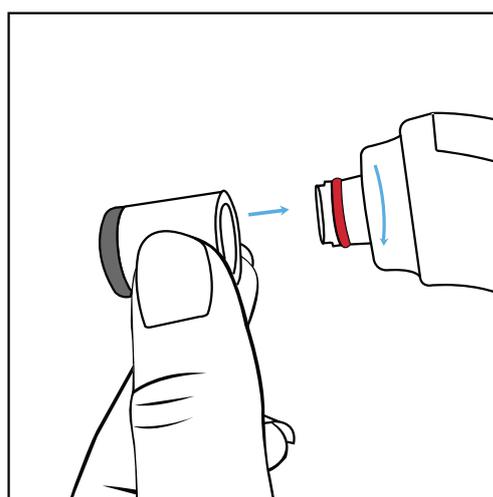

(Fig 7)

6. Turn on the device by pressing the power button once. The vacuum level LED will light red indicating GLO-ME is set on High.
7. Adjust the vacuum level by pressing on the power button to switch between:
 - a. Red – High
 - b. Orange – Medium
 - c. Green – Low
 Your first GLO-ME treatment should be set on LOW.
8. Before use, test on a small and concealed skin area and check if any skin irritation appears.
9. At the end of the procedure remove the filter and dispose of it.
10. Clean the tip according to cleaning instructions.
11. For post treatment care it is recommended to apply sunblock.

TREATMENT TECHNIQUES

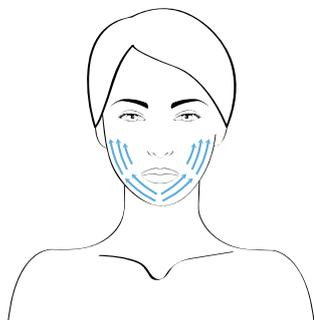

(Fig 8)


Lower Face

Start treatment at the chin, move along the jaw and then towards the ear using gentle upward strokes. Move the tip higher to treat the whole cheek area (Fig 8).

Move to the other side of the face and repeat as above (Fig 8).

Forehead

Treat the forehead starting at its center and use gentle outward strokes on one side and when that area is covered, move to the other side of the forehead (Fig 9).

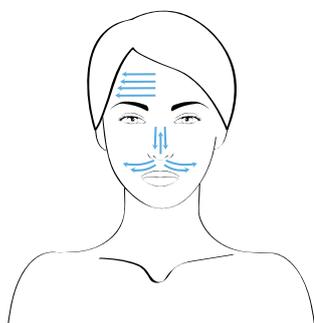

(Fig 9)

Nose

Treat the nose starting on one side and gently gliding the GLO up and down. Then switch to the other side (Fig 9).

Upper Lip Area

Treat the area above the upper lip starting from the center using gentle sideway strokes on one side and then move to the other side (Fig 9).

Notes:

- When taking the tip off your skin, do it gently and gradually so it will not irritate the skin.
- **Do not use directly above or below the eyelids and never on the eye socket.**
- Do not overuse as it may injure and harm your skin.
- If you experience any discomfort, discontinue use immediately and consult a physician.
- You can also treat the neck, décolleté or hands using the same type of outward strokes technique making sure to cover all areas.

TREATMENT PROTOCOL

- Use every other day. Make sure skin is completely healed before applying the next treatment.
- During treatment, do not move the GLO-ME with back and forth motions (move in one direction).
- Do not move over the same area more than 3 times.
- Stopping or moving slowly may cause scratching or bruising of the skin.
- Do not use excessive pressure.
- Excessive use may cause pigmentation or dermatitis depending on the individual's skin condition and characteristics. If you observe any skin problems after use, stop the treatment immediately and consult with a dermatologist.



CLEANING & MAINTENANCE

Shut off the device before cleaning or maintenance.

Cleaning

Tip: Soak the diamond peel tip in hot water (50~70°C) for approximately 5-10 minutes, and rinse with running water. Use a brush with soap (or cleanser) to remove any remaining dirt. Dry thoroughly.

Do not place the tip in alcohol or boiling water (100°C).

Device: Clean with a damp cloth. Do not use soap or alcohol.

Storage & maintenance after use

- Avoid storing and using in high humidity or wet areas.
- Do not immerse in any liquids.

CUSTOMER SUPPORT & WARRANTY

Warranty Terms

This product has been subjected to strict quality control performed by the manufacturer. For this reason, the GLO-ME™ distributor in each country guarantees proper function of your product for 12 months from date of actual purchase. To use the device properly, follow all instructions outlined in the user manual and you will enjoy the product to your ultimate satisfaction.

This warranty ensures that the distributor will provide support for the GLO-ME™ device, in case of failure due to manufacturing defects. In such case, you shall be entitled to make full use of the benefits of this warranty. This warranty is valid for the GLO-ME™ device marketed by the distributor in your country.

The warranty is not valid when:

1. Amendments in the sales receipt or invoice, or the absence of any of these original documents.
2. Damage or deterioration of the device due to improper use or if the device has been modified or repaired by an unauthorized third party.
3. Damage of the product due to defects incurred during transportation, such as bumps, bangs or improper handling.
4. Damage due to humidity or any kind of liquids, foreign bodies inside the device or the use of unapproved, unadvised accessories with the device contradictory to the user manual.
5. When the device is not used in accordance with the respective user manual or with its supplied accessories.

6. Partial or total damage or loss of the device is due to natural disasters (such as earthquakes, floods, fire, lightning, etc.), force majeure or accident.
7. Damage to the device after being subjected to high dust environments, humidity or application of excessive voltage from the power inlet.
8. If the product's serial number has been amended or erased.
9. To cosmetic damage, including but not limited to scratches and dents.

GL051400001



Device S/N location

TECHNICAL SPECIFICATIONS

Dimensions (mm):	44(W) x 70(D) x 170(H)
Weight:	100g (including batteries)
Power supply:	1.5V, AAA size battery x 2
Max noise level:	72 dB

Disclaimer: Users of this device use it at their own risk. Neither POLLOGEN BEAUTY nor its retailers assume any responsibility or liability for any injuries or damages, physical or otherwise, resulting, directly or indirectly, from the use of this device.

Model may be changed for improvements without notice.